## Instructions for use Ergoject Top-Loaded

EN	Instructions for use	Page	2
DE	Gebrauchsanweisung	Seite	5
AR	تعليمات الاستخدام	صفحة	8
BG	Инструкции за употреба	страница на книгата	11
CS	Pokyny k použití	S stránka knihy	14
DA	Brugsanvisning	Bog side	17
EL	Οδηγίες χρήσης	σελίδα βιβλίου	20
ES	Instrucciones de uso	Página	23
ET	Kasutusjuhend	Raamatu leht	26
FI	Käyttöohjeet	Kirjan sivu	29
FR	Mode d'emploi	Page	32
HR	Upute za upotrebu	Stranica knjige	35
HU	Használati utasítás	Könyv oldala	38
IT	Istruzioni per l'uso	Pagina	41
KO	사용 지침	예약 페이지	44
LT	Naudojimo instrukcija	Knygos puslapis	47
LV	Lietošanas instrukcija	Grāmatas lapa	50
NL	Gebruiksaanwijzing	Pagina	53
NO	Bruksanvisning	Bokside	56
PL	Instrukcja użytkowania	Strona książki	59
PT	Instruções de utilização	Página	62
RO	Instrucţiuni de utilizare	Pagina	65
RU	Инструкция по применению	Страница книги	68
SK	Návod na použitie	Stránka knihy	71
SL	Navodila za uporabo	Stran knjige	74
SR	Uputstvo za upotrebu	Stranica knjige	77
SV	Bruksanvisning	Sida	80
TR	Kullanma talimatları	Kitap sayfası	83
UK	Інструкція з експлуатації	Сторінка книги	86

### 1. Scope of application

These instructions for use apply to the following medical devices:

Reference number	Product name
ERT1600	Ergoject 1.6-TL Injector
ERT1800	Ergoject 1.8-TL Injector
ERT2000	Ergoject 2.0-TL Injector
ERT2100	Ergoject 2.1-TL Injector
ERT2200	Ergoject 2.2-TL Injector
ERT2400	Ergoject 2.4-TL Injector
ERT2600	Ergoject 2.6-TL Injector
ERT3000	Ergoject 3.0-TL Injector

### 2. Scope

Intended use: Folding and injecting of an intraocular lens (IOL) during eye surgery.

Indication: Visual problems that require implantation of IOL.

Contraindications: There are no specific contraindications for the injectors, apart from general contraindications related to ophthalmic surgery.

  User Group: The medical device must only be handled by healthcare professionals and applied to the patient by the surgeon. Use by unqualified personnel could result in patient injury.

 Patient Group: A patient-related, careful preoperative assessment, with well-founded clinical judgment and risk-benefit assessment, is the responsibility of the surgeon.

Clinical benefit: IOL injector allows IOL implantation through small incisions.

### 3. Notes on safety

The suitability of the intraocular lens (IOL) in combination with the medical device (IOL injector) must be tested and approved in advance by the intraocular lens manufacturer.

  **Warning**  
Using a damaged or non-sterile medical device can lead to infection / endophthalmitis in the patient.  
Do not use the medical device if there are signs of damage to the medical device itself or to the sterile packaging.

 **Warning**  
The medical device must not be used after the expiry date. Use after the expiration date may result in infection / endophthalmitis of the patient.

 **Warning**  
Touching the endothelium with the silicone cushion or IOL can damage the corneal endothelium.

 **Warning**  
Appropriate surgical technique is the responsibility of the individual surgeon. The surgeon has to assess the suitability of the respective procedure on the basis of his / her training and experience.

 **Warning**  
The medical device is intended for single use only. Reuse or inappropriate reprocessing may lead to serious adverse effects on the health and safety of the patient.

 - Reuse may result in infection / endophthalmitis, TASS or allergic reaction due to cleaning residues.

 - The performance of the medical device may be compromised. Heat or the cleaning process will destroy performance properties (e.g. mechanical properties and gliding properties) of the medical device.

 **Caution**  
In order to avoid capsule rupture, the IOL must be injected slowly and carefully into the eye.


#### 4. Reporting to manufacturer and authorities

Serious incidents occurring in connection with the device must be reported to the manufacturer and the competent authority in which the user and/or patient is established.

#### 5. Application instructions for the single-use IOL injector

1. Open the blister in a sterile environment and remove the sterile injector (1).
2. Hold the injector so that the rear wing of the loading chamber can be guided with the index finger of your left hand (2).
3. Fill the cartridge tip and the loading chamber sufficiently with ophthalmic viscosurgical devices (OVD) (3).  
  
It is recommended to flush the cartridge tip and the loading chamber with Balanced Salt Solution (BSS) before the application of OVD.  
  
Allow the OVD to react (30s). OVD can lose their lubricating properties if they are in contact with air for too long. Therefore, the IOL should be injected without delay after loading.  
  
Hydrophobic IOL types: The application of OVD is required. The exclusive use of BSS is not permitted. In addition, apply a drop of OVD to the silicone cushion (3).
4. Open the loading chamber so that the IOL can be easily loaded. Position the IOL in the middle of the loading chamber with a concave orientation (4).
5. Use the sterile forceps to apply slight pressure to the IOL optics (5), this ensures that the IOL folds in a concave direction. For IOLs with C-loop haptics, align the haptics with sterile forceps in the direction of the optics.
6. At the same time, close the wings of the loading chamber until the „Click-Lock“ mechanism engages (6)..
7. Push the injector plunger forward.  
  
If the silicone cushion is inclined and cannot be inserted into the closed loading chamber, the silicone cushion can be aligned with the forceps.  
  
Pull back the injector plunger a few millimetres and then push it forward again. This step ensures that the IOL is always gripped correctly.
8. Push the injector plunger forward until the rear pusher plate is flush with the injector housing or until the wheel moves (7).  
  
The IOL is now loaded and ready for injection.
9. Use the small lever on the injector housing to choose whether the injector should be operated in forward and reverse mode (lever in middle position) or in forward mode only (lever to the left) (8). In the forward mode only, the return movement of the injector plunger is blocked. This is advantageous if the IOL and silicone cushion generate a strong recoil, which makes it difficult to operate the injector wheel (especially with models for small incisions). The activation of the exclusive advance mode is signalled by a slight click when the injector plunger is pushed forward. At any time during the injection process, it is possible to switch back and forth between forward and reverse mode and exclusive forward mode by simply adjusting the lever.
10. With your index finger, slowly pull the wheel back to advance the IOL (9).  
  
If the rear haptic is wedged between the injector plunger and the cartridge tip after the plunger has been fully advanced, set the injector to forward and reverse mode (lever in the middle position) and pull the injector plunger back a few millimetres until the haptic is exposed. Activate exclusive forward mode again if necessary and inject the rear haptics in a second push.
11. Slowly inject the IOL into the eye.
12. The correct exit of the IOL can be supported by turning the injector slightly. Only push the injector plunger forward until the IOL has emerged completely.
13. If necessary, assist the IOL with the help of a suitable positioning hook during the exit and bring the IOL into its final position.
14. Thoroughly remove the viscoelastic material from the eye and IOL using standard irrigation and aspiration techniques.
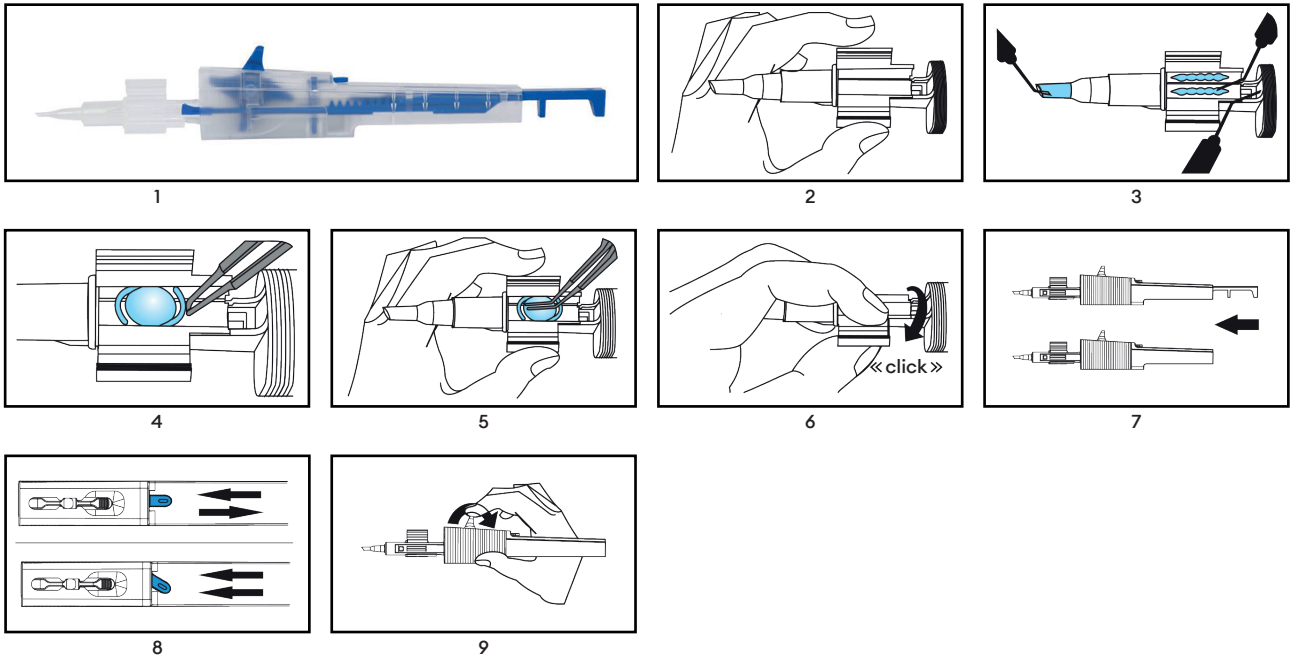
## 6. Disposal

### Warning

The medical device may be contaminated with potentially infectious substances of human origin after use.

Dispose of the medical device and packaging after use according to the applicable guidelines of biological hazardous waste.

## 7. Symbols

	Sterilized using ethylene oxide		Keep dry
	Number of medical device in packaging unit		Do not use if package is damaged and consult instructions for use
	Single sterile barrier system		Keep away from sunlight
	Single sterile barrier system with protecting packaging outside		Do not re-sterilize
	Date of manufacture YYYY-MM-DD and country of manufacture		CE symbol with number of the notified body
	Do not re-use		Authorized representative in the European Community / European Union
	Reference number		Open here
	Lot number		Unique device identifier
	Expiry date		Medical device
	Manufacturer		Consult electronic instructions for us <a href="http://www.medicel.com/ifu">www.medicel.com/ifu</a>
	To be used by qualified personnel only		Caution
	Recycling		

### 1. Scope of application

These instructions for use apply to the following medical devices:

Referenznummer	Names des Produkts
ERT1600	Ergoject 1.6-TL Injector
ERT1800	Ergoject 1.8-TL Injector
ERT2000	Ergoject 2.0-TL Injector
ERT2100	Ergoject 2.1-TL Injector
ERT2200	Ergoject 2.2-TL Injector
ERT2400	Ergoject 2.4-TL Injector
ERT2600	Ergoject 2.6-TL Injector
ERT3000	Ergoject 3.0-TL Injector

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